



426 ADIMAL™ Digital Thermometer

Our newest digital thermometer with a fresh face.

The Adimals™ 426 Thermometer features

- Jumbo display
- Range 89.6°F-109.4°F $\pm 0.3^\circ\text{F}$ or $32^\circ\text{C} - 43^\circ\text{C} \pm 0.2^\circ\text{C}$ depending upon scale selection**
- Replaceable 1.55v (LR41) type battery provides up to 4,000 measurements
- Clear plastic scabbard carry case
- Water resistant

Ordering Information

Item	Description
426*	Adimal™ Digital Thermometer / Dual Scale *Specify Bear-B, Frog-F, Koala-K, or Panda-P




2150 ADIMAL™ Digital Fingertip Pulse Oximeter

Our newest pediatric fingertip pulse oximeter with a fresh face!

The 2150 features

- Reads SpO_2 and Pulse Rate
- Accuracy: $\text{SpO}_2 \pm 3\%$ / HR ± 2 BPM
- Six changeable display modes
- 10 display brightness settings
- Child-friendly instrument face
- Fits fingers with a thickness of 0.2" to 0.63" (5.0mm to 16.0mm)
- Carry lanyard
- Soft carry case
- Two "AAA" batteries
- Instruction booklet
- Two-year limited warranty

6 changeable display modes

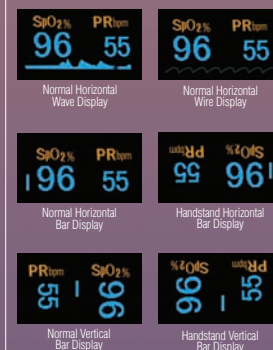

Carry Case and Lanyard



Ordering Information

Item	Description
2150	Adimal™ Fingertip Pulse Oximeter
2100-1	Pulse Oximeter Carry Lanyard





Adscope™ 618 ADIMAL™ Scope

Our newest pediatric multifrequency scope with a fresh face!

The 618 series features

- Distinctive, molded-resin animal face chestpiece optimally sized (1½") for pediatric use
- Multifrequency design switches between bell-like low frequency and diaphragm-like high frequency performance simply by altering pressure
- Non-chill bell and diaphragm retaining rim for patient comfort
- Clinical headset boasts reinforcing yoke molded into the flexible 22" PVC tubing
- Lightweight aluminum binaurals in matching satin finish
- Adsoft™ threaded PVC eartips for the ultimate in wearing comfort and acoustic seal
- Seven interchangeable animal face snap-ons included:
 - Panda
 - Frog
 - Koala
 - Deer
 - Monkey
 - Tiger
 - Bear
- Weighs: 3 oz. (85 g.)
- Overall length 30"
- Five-year warranty



Multifrequency Diaphragm

Ordering Information

Item	Description
618*	Adscope™ 618 Adimals™, Pediatric (Specify Color)

Color Choices



NOTE: To ensure correct color when ordering.
Add Color suffix to model number (i.e., 618**SF**)



www.adctoday.com

